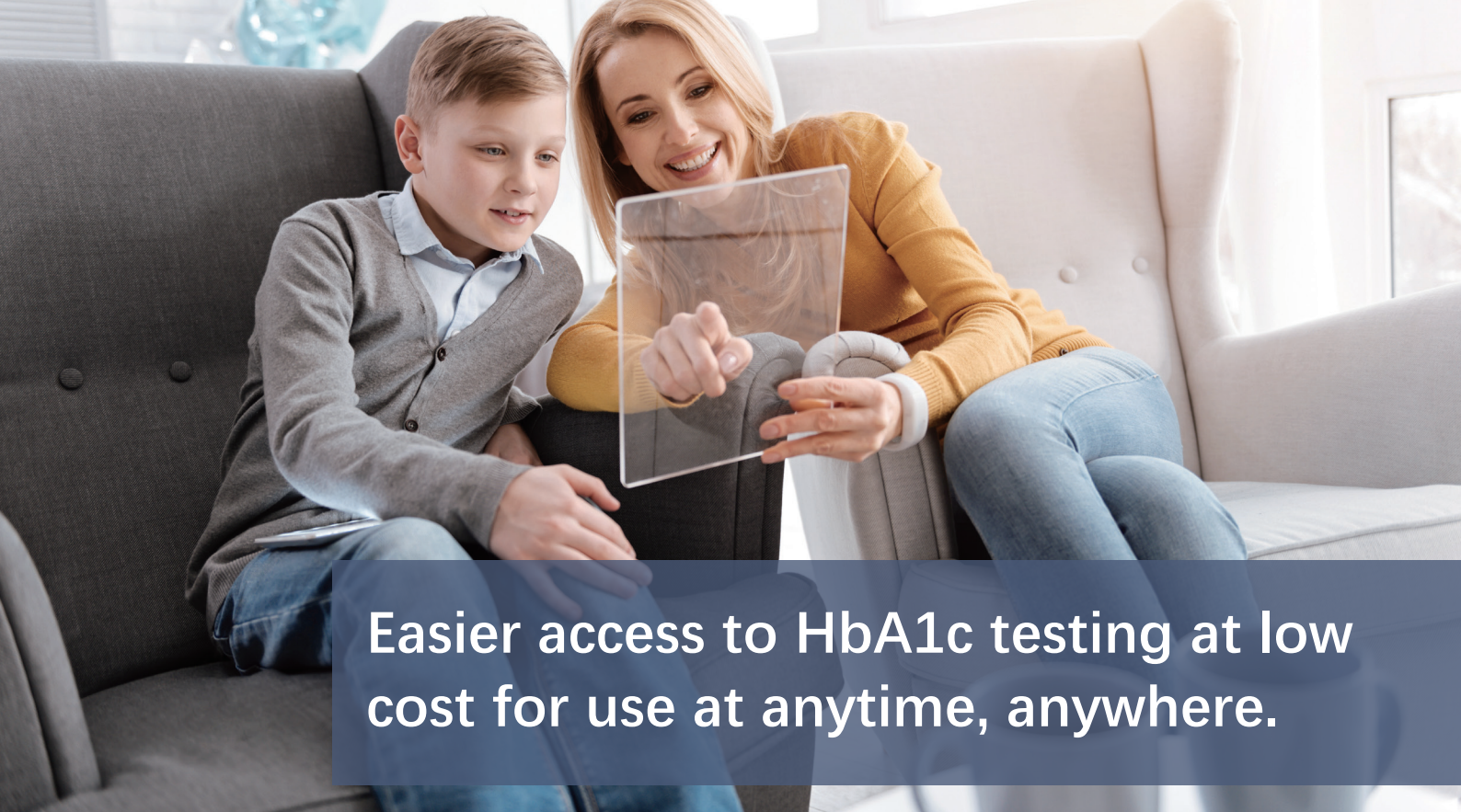




Easier access to HbA1c testing at low cost for use at anytime, anywhere.

GluCoA1c Blood Glucose and Glycohemoglobin Analyzer

HbA1c and blood glucose testing for professional and self-use

Enhanced testing function

Blood Glucose & HbA1c Testing
-2 parameters in 1 device.

Easy and simple to handle

Handheld HbA1c device;
Voice prompt operation, easy to perform.

Fast & Reliable testing result

Two results within 5 minutes; comprehensive patient diabetic management at any site.



Product Specifications

Parameter	Specification
Principle	Boronate affinity chromatography
Test Items	HbA1c/Blood glucose
Sample Type	Capillary/Venous blood for HbA1c test Capillary blood for blood glucose test
Sample Volume	Around 3 μ L
Testing Time	5 minutes
Metric	NGSP (%) /IFCC (mmol/mol)/eAG (mg/dL or mmol/L)
Test Range	4.0%-14.0%
CV	< 3%
Strip and Reagents Storage Condition	18 months at room temperature
Analyzer Size	126 x 64.6 x 23.5 mm
Analyzer Weight	99.5 g
Analyzer Power Supply	3.7V Lithium battery
Working Condition	Temperature: 10-40 $^{\circ}$ C ; Humidity: 10%-90%
Printer	Support to connect the external printer
Data Storage	1000
Date Transfer and Connectivity	USB/Bluetooth

Order Information

Product Name	Catalog No.	Components
GluCoA1c Blood Glucose and Glycohemoglobin Analyzer	GCA-M11 GCA-M31 GCA-M41	1 GluCoA1c Analyzer 1 Battery 1 Lancing Device 1 User Manual 1 Quick Start Guide 1 Cleaning Guide 1 Warranty Card

Product Name	Catalog No.	Components
Glycohemoglobin Test Kit	GCA-S12 GCA-S22 GCA-S32	12 HbA1c Test Strips 12 Samplers 1 Buffer A 1 Buffer B 1 Code Chips 1 Package Insert

Product Name	Catalog No.	Components
GluCoA1c Blood Glucose Test Strips	GLU-S12 GLU-S22 GLU-S32	50 Test Strips 1 Package Insert



GluCoA1c
Blood Glucose
and Glycohemoglobin Analyzer



BioHermes® BioHermes.com

Wuxi BioHermes Bio & Medical Technology Co., Ltd.
Tel: + (86) 510-85385897 Ext: 2091
Fax: +(86) 510-85387563
E-mail: info@biohermes.com.cn